Total Remove Kit




Most easy way to remove broken screw & fixture by using Total Remove Kit.



HNTRK-300


For detail of information and video by scanning QR code.



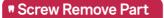

Implant Guide Using screw remove drill/tap/bar depend on situation in reverse side to remove broken screw with guide.



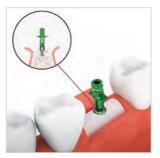





Fixture Remove Part

Fixture Remove



Using fixture remove take fixture and then reverse side to remove fixture by attached ratchet wrench and adapter.



■ Cover Screw Remove



Using cover screw remover when you find difficulty to open the closing screw.



■ Instruments

Adapter



12 HHSA40

Hand tool to use fixture remove for handy type.



Square Wrench





Hand tool to use fixture remove for high torque.

